

PressureGuard® APM².0

PressureGuard® APM² with digital control unit



Now made with
LifeSPAN™
fortified healthcare fabric
bleach-resistant fabric cover
(see page 15 for details)

Therapy Provided

- ✓ Powered Pressure Redistribution
- ✓ Alternating Pressure
- ✓ Lateral Rotation

The APM² has been an industry-leading air therapy solution for more than a decade, based in large part on its documented superiority vs. typical, roll-up treatment surfaces. Newly upgraded cover design and control unit provide even greater immersion, envelopment, shear reduction, treatment flexibility, caregiver ergonomics and user comfort.

QUICK FACTS

Indications:

- Treatment of Stage 1-4 pressure injuries
- Patients that can't/won't reposition frequently
- Can be positioned off pressure injury in at least two positions¹
- Prevention for high risk patients

¹ In situations where positioning on the pressure ulcer cannot be avoided (e.g. the individual has multiple ulcers on multiple surfaces), limit the amount of time the individual is positioned on the ulcer. Source: NPAUP/EPUAP Pressure Ulcer Prevention and Treatment, Clinical Practice Guideline, p 66-67

– New, Upgraded Features –

Digital Multi-function Control Unit



Lightweight and quiet, it delivers an unmatched combination of treatment flexibility, ease of use, and dependability.

- Low Pressure indicator light and selectable audible alarm.
- Power Failure indicator light and selectable audible alarm.
- Simple, indexed comfort adjustment.
- Timed (20-minute) "Auto Firm" mode for transfers, ADLs, & CPR.
- "Float" mode provides powered flotation therapy.
- "Alternate" mode changes loading in 10-minute cycle; controlled at mattress for choice of lateral rotation or alternating pressure therapy.




Shear Transfer Zones® Cover Design (US Patent Number 8,438,682)

In addition to the Geo-Matt® segmented top surface, Span-America's exclusive Shear Transfer Zones® provides an additional measure of shearing protection in the form of the silicone-coated, shear-minimizing fabric bands located on the underside of the bi-directional stretch cover.

This Shear Transfer Zones design helps:

- Prevent heels, sacrum, and scapula from "digging into" surface.
- "Glide" user back to original position following HOB elevation.
- Protect against damaging effects of micro shear, macro shear, and rotational (pivot-induced) shear.
- Provide patient stability by transferring shear to more shear-tolerant anchor points.

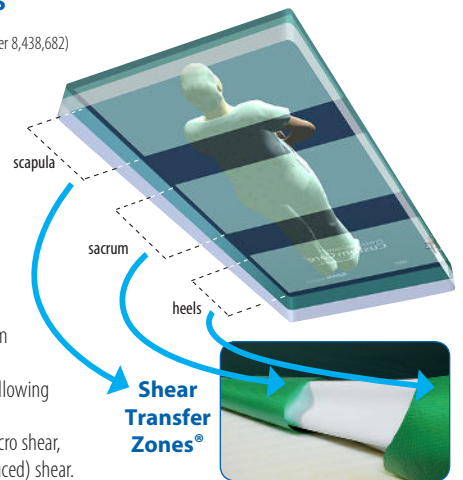

Item #	Description	Dimensions
5875LR-29	APM ² with digital control unit	75"L x 35"W x 7"H
5880LR-29	APM ² with digital control unit	80"L x 35"W x 7"H
5884LR-29	APM ² with digital control unit	84"L x 35"W x 7"H
* AP28039-29	APM ² with digital control unit	80"L x 39"W x 7"H
* AP28439-29	APM ² with digital control unit	84"L x 39"W x 7"H
* 3715629	APM ² with digital control unit	80"L x 42"W x 7"H
5900	Digital control unit only	

* Allow up to 10 working days lead time for this size.

Control Unit: Weight: 6.3 lbs. Dimensions: 12" x 7.5" x 5"
Voltage: 120 AC Max current: 1.0 Amp
Frequency: 60 Hertz UL Listing: Certified to IEC 60601-1-1

Medicare Approved: E0277 (Group II)
Weight limit: 350 lbs
Mattress weight: approx. 20 lbs.
Warranty: 18 months (all components)

All models are certified for conformance to NFPA 101 (Life Safety Code) – ASTM E1590; Cal Tech Bulletin #117 and 16 CFR Parts 1632 and 1633 by an independent testing organization.

WARNING: Cancer and Reproductive Harm – www.P65Warnings.ca.gov



SURFACES