



CARLO Classic Floor Lift

Safe and Effortless Patient Transfers



Providing a High Level of **Safety and Comfort**

Characterized by a very light and stable aluminum construction, CARLO Classic supports caregiver and resident safety in daily lifting, transfer and repositioning tasks.

Featuring a 6-point carry bar that can be easily lowered to the floor for lifting residents in slip and fall accidents, the low height of the carry bar can also facilitate transfers to bathtubs and beds. Featuring an Anti-Rotation Mechanism (ARM), CARLO Classic removes any involuntary rotation of the carry bar during transfers and reduces uncomfortable swivel movements, providing the resident an extra feeling of safety.

CARLO Classic is available with a safe working load of 230 kg (507 lb) for bariatric residents with low mobility. Also available with low legs to accommodate lower beds, stretchers and furniture.

Battery and charger are backwards compatible across all Handicare by Beka mobile lifts.

Caregiver and Resident Benefits




Carrybar with anti-rotation mechanism (ARM) and a free rotation area of 360°.



Cable remote control with display of the charging condition of the battery, failure diagnostics, service.



Four high-quality double castors, two of them with brakes, provide easy movement.

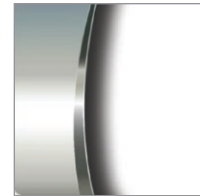

Light construction. Chassis, lifting column and components made of 100% aluminum.



Electrically spreadable legs with low height to pass underneath baths and beds.



Multi-purpose 6-point carrybar allows an easier and comfortable transfer of the client.



Electrical height adjustment with lowering of the carrybar to the floor. Also includes manual emergency lowering system.



Wide choice of loop slings.

CARLO Classic Floor Lift

Specifications

Safe Working Load (SWL)

CARLO Classic: 185 kg (408 lb)
CARLO Classic 230: 230 kg (507 lb)

Dimensions (L x W x H)

CARLO Classic

Carrier frame, length: 1290 mm (50.8 in)
Carrier frame, height: 125 mm (4.9 in)

Inner width, closed: 540 mm (21.3 in)
Inner width, spread: 1040 mm (40.9 in)

Outer width, closed: 700 mm (27.6 in)
Outer width, spread: 1180 mm (46.5 in)

Total width of lift, carrier frame closed:
700 mm (27.6 in)

Height: 1340 – 2170 mm (52.8 – 85.4 in)

CARLO Classic 230

Carrier frame, length: 1400 mm (55.1 in)
Carrier frame, height: 125 mm (4.9 in)

Inner width, closed: 540 mm (21.3 in)
Inner width, spread: 1090 mm (42.9 in)

Outer width, closed: 700 mm (27.6 in)
Outer width, spread: 1230 mm (48.4 in)

Total width of lift, carrier frame closed:
700 mm (27.6 in)

Height: 1340 – 2170 mm (52.8 – 85.4 in)

For additional technical data, please refer to the user manual.

Lifting Range

1400 mm (55.1 in)

Castors

100 mm (3.9 in)

Power Supply

Removable 24V battery

Product Weight Without Packaging

55 – 65 kg (121 – 143 lb)

Part Numbers

924003208: CARLO Classic
On board charging and dual controls

924002508: CARLO Classic with Low Legs
On board charging and dual controls

924003218: CARLO Classic
Includes two batteries, wall mount charger and dual controls

924003008: CARLO Classic 230
On board charging and dual controls

924003018: CARLO Classic 230
Includes two batteries, wall mount charger and dual controls

Accessories / Options

CARLO Classic

Description	Part No.
Conversion kit for retrofitting a patient scale	924003400
Calibrated patient scale	924003300

CARLO Classic 230

Description	Part No.
Conversion kit for retrofitting a patient scale	924001800
Calibrated patient scale up to 250 kg (551 lb)	925001208

For all floor lifts

Description	Part No.
Spare battery, 24V, blue	273225
Spare battery, 24V, orange	921070602

Are you ready to make safe and effortless resident transfers part of your everyday?

To learn more, call 1-888-637-8155 or visit www.handicare.ca / handicareusa.com

HANDICARE

485 Millway Avenue, Unit 2, Concord, ON, L4K 3V4

Toll Free: 877.304.5438 | Fax: 416.260.0256

customerservice.canada@handicare.com

www.handicare.ca

10888 Metro Court, St. Louis, MO 63043

Toll Free: 866.891.6502 | Fax: 314.692.7858

customerservice.northamerica@handicare.com

www.handicareusa.com