



NORA Sit-to-Stand Lift

Light Flexible Standing Aid for Daily Care



Providing a **High Quality of Life** for Your Residents

NORA is designed to assist residents who can stand for short periods of time but require a safe and comfortable transfer to and from bathing or toileting. For shorter residents NORA offers additional attachment points for slings so residents can stand upright independent of their height, the flexible knee pads can also be adjusted to suit height requirements. Residents can also boost circulation and rehabilitate with NORA; the targeted safe raising system allows residents to slowly put their full vertical load onto their legs improving musculoskeletal systems.

All electrical components are splash water resistance.

Battery and charger are backwards compatible across all Handicare by Beka mobile lifts.

Caregiver and Resident Benefits




The ergonomic handles allow for an optimum holding position for steering.



Two sling attachment points for different patient heights.



Cable remote control with display of the charging condition of the battery, failure diagnostics, service.

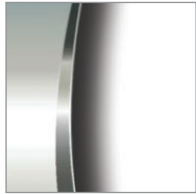

Light construction. Carrier frame, lifting column and add-on components made of 100% aluminum.



Adjustable, extra soft double knee pads.



Four high-quality double castors, two of them with brakes, provide easy movement.



Electrically spreadable carrier frame.



Can also be used for showering the seat area.



Specially designed, high-quality lifting slings available in different sizes.

NORA Sit-to-Stand Lift

Specifications

Safe Working Load (SWL)

200 kg (440 lb)

Dimensions (L x W x H)

Length: 1027 mm (40.4 in)

External width, closed carrier frame : 700 mm (27.6 in)

Internal width, spread carrier frame: 897 mm (35.3 in)

Height: 1168 mm (46.0 in)

Height, carrier frame: 126 mm (5.0 in)

Lifting Range

Minimum: 898 mm (35.0 in)

Maximum: 1726 mm (68.0 in)

Castors

Front: 100 mm (3.9 in)

Rear: 125 mm (4.9 in)

Power Supply

Removable 24V battery

Product Weight Without Packaging

51 kg (112 lb)

For additional technical data, please refer to the user manual.

Part Numbers

921071008: NORA Sit-to-Stand Lift

On board charging and dual controls

Accessories / Options

NORA Stand-up / Raising Sling

Description	Part No.
NORA Sit to Stand Sling – Small	921070053
NORA Sit to Stand Sling – Medium	921070103
NORA Sit to Stand Sling – Large	921070203
NORA Sit to Stand Sling – Extra Large	921070303

NORA Double Knee Pad

Description	Part No.
Extra soft double knee pad for NORA	921070408

NORA Carry Bar Loop Pin Kit

Description	Part No.
NORA Carry Bar loop pin kit (pair)	921071030

Spare Battery

Description	Part No.
Replacement battery 24V (orange)	921070602

Are you ready to make safe
and effortless resident transfers
part of your everyday?

To learn more, call 1-888-637-8155 or
visit www.handicare.ca / handicareusa.com

HANDICARE

485 Millway Avenue, Unit 2, Concord, ON, L4K 3V4

Toll Free: 877.304.5438 | Fax: 416.260.0256

customerservice.canada@handicare.com

www.handicare.ca

10888 Metro Court, St. Louis, MO 63043

Toll Free: 866.891.6502 | Fax: 314.692.7858

customerservice.northamerica@handicare.com

www.handicareusa.com