



SENTA PUR L Bath and Shower Lift

Innovative Combination Lift for Bathing and Showering



Easily Converts Into a Stretcher or Bath Lift

The SENTA PUR L is an innovative bath lift that quickly converts into a stretcher style lift with the addition of optional “wing” modules, offering unparalleled leg support to residents with limited mobility. The SENTA PUR L can also lift residents directly from beds. The deep carrier frame allows for easy pass through for simplified pick-up onto the lift. The unique combination style design of the SENTA PUR L eliminates the need for purchasing an additional lift, caregivers can easily switch from seated position to the lying position providing your resident with versatility in care. Fits any type of existing care tub, the SENTA PUR L is extremely easy to maneuver and can be operated by a single caregiver.

Battery and charger are backwards compatible across all Handicare by Beka mobile lifts.

Caregiver and Resident Benefits



Handheld operation with problem indication and battery level display.



Digital patient scale. (optional)



Seat area with large, lockable hygiene opening.



Easy battery access and removal.




The PUR padding is easily slid onto the stainless steel frame and can be removed at any time.



Four high-quality double castors, two of them with brakes, provide easy movement.



The carrier frame is made from high-quality steel and formed for self-entry.

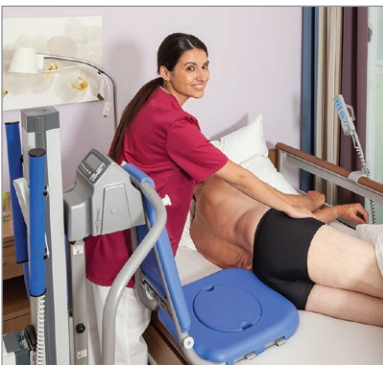

Easily lift residents directly from beds with easy pass through.



Side rail for increased safety plugs in easily and features safety lock.



Two removable stretcher lift wings comfortably convert the seat lift into a full stretcher lift. The two modules are easily clicked into an existing safety adapter at the seat tray. Also available with a foot rest. (optional)

SENTA PUR L Bath and Shower Lift

Specifications

Safe Working Load (SWL)

150 kg (331 lb)

Dimensions (L x W x H)

Length: 860 mm (33.9 in)

Width: 718 mm (28.3 in)

Height: 1375 mm (54.1 in)

Height, seat: 475 mm (18.7 in), minimum

Height, chassis: 125 mm (4.9 in)

Free space under carrier frame: 25 (1.0 in)

Lifting Range

475 – 1175 mm (18.7 – 46.3 in)

Castors

Two double castors, front: 100 mm dia. (3.9 in)

Two double castors, back: 125 mm dia. (4.9 in)

Power Supply

Removable 24V battery

Protection Class

IPX4

Product Weight Without Packaging

63 kg (139 lb)

For additional technical data, please refer to the user manual.

Part Numbers

920101008: SENTA PUR L Bath and Shower Lift
Dual control box with onboard charging

Accessories / Options

Description	Part No.
Digital patient scale	920104008
Stretcher wings (2 pieces)	920103008
Side rail	920103108
Foot rest	920103115



HANDICARE

485 Millway Avenue, Unit 2, Concord, ON, L4K 3V4

Toll Free: 877.304.5438 | Fax: 416.260.0256

customerservice.canada@handicare.com

www.handicare.ca

10888 Metro Court, St. Louis, MO 63043

Toll Free: 866.891.6502 | Fax: 314.692.7858

customerservice.northamerica@handicare.com

www.handicareusa.com

Are you ready to invest in innovation
and simplify bathing for your residents?

To learn more, call 1-888-637-8155 or
visit www.handicare.ca / handicareusa.com