



# NORA Pro Sit-to-Stand Lift

Next Generation Design and Lifting Comfort



## Designed for a More **Natural, Ergonomic, Standing-Up Motion** with Excellent Comfort

NORA Pro sets new standards in both comfort and safety in sit-to-stands. The centerpiece of the newly designed lift is the unique 3D-Comfort-System consisting of a large, flexible, memory-foam cushion that offers excellent leg support that surrounds the client in comfort throughout the transfer. Research guided the R&D teams to feature a dementia-sensitive colour design supporting clients with orientation difficulties while also providing them highly visible contact points.

Compact and lightweight, NORA Pro features an aluminum construction with excellent maneuverability, making it suitable for both toileting and bathing in limited spaces. The additional sling attachment points also ensure a comfortable transfer for clients with a body height under 5.3 feet.

All electrical components are splash water resistance and rated IPX4.

# Caregiver and Client Benefits



Ergonomically shaped handles for a secure hold; additional user grips underneath bar.



Features four sling attachment points for tall and short clients. The additional attachment points also prevent unpleasant sling pressure under the arms.



Equipped with an LCD display showing battery status, cycle counter, preventative maintenance indicator, and SWL alerts. The unit features lift controls in addition to the hand control.




The optional Leg Separator provides additional stability and can support the care staff in changing incontinence briefs.



The fixed safety belts of NORA Pro are equipped with a retractable mechanism for added safety.



The horizontal footplate is easily detachable for easy cleaning.

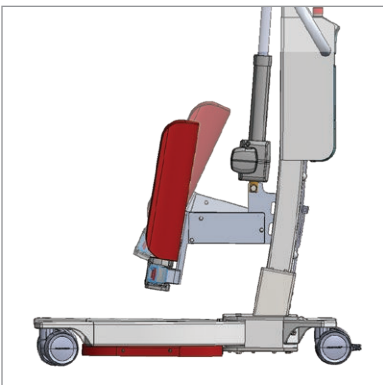

The 3D-Comfort-System enables ergonomic lifting through its two swivel axles, adapting to the natural sit-to-stand movement patterns of the legs to provide maximum comfort.



Leg support, grip handles and the footplate of NORA Pro are coloured in a pleasant red and provides high perceptual value for your clients. Red is most likely to be noticed by people suffering from dementia.\*



The multi-layer memory foam padding offers a secure grip in the cushion, the user's legs are comfortably surrounded by the memory-foam during the lifting process. The large cushion supports different legs sizes and shapes.

\*Studies show that red is best perceived by people suffering from dementia and with gradually declining eyesight. NORA Pro therefore supports clients with its high contrast colour concept.

# NORA Pro Sit-to-Stand Lift

## Specifications

### Safe Working Load (SWL)

205 kg (451 lb)

### Dimensions (L x W x H)

Length: 944 mm (37.17 in)

Width, minimum : 704 mm (27.72 in)

Inside width legs, chassis spread: 910 mm (35.82 in)

Height, minimum: 1091 mm (42.95 in)

Height, maximum: 1608 mm (63.31 in)

Height, carrier frame: 1091 mm (42.95 in)

### Lifting Range

Minimum: 898 mm (35.35 in)

Maximum: 1614 mm (63.54 in)

### Castors

Front: 100 mm (3.93 in)

Rear: 100 mm (3.93 in)

### Leg Clearance

129 mm (5.08 in)

### Power Supply

Removable 24v Lithium-Ion battery

### Product Weight Without Battery

49 kg (109 lb)

### References

\*Wijk, H., Berg, S., Sivik, L., & Steen, B. (1999). Colour discrimination, colour naming and colour preferences among individuals with Alzheimer's disease. *International Journal of Geriatric Psychiatry*, 14(12), 1000–1005.

\*Dunne, T.E., Neargarder, S.A., Cipolloni, P. B., & Cronin-Golomb, A. (2004). Visual contrast enhances food and liquid intake in advanced Alzheimer's disease. *Clinical nutrition (Edinburgh, Scotland)*, 23(4), 533–538.

## Part Numbers



**921078000:** NORA Pro RUBIN (Clip Connection)  
Dementia-sensitive colour concept

**921078001:** NORA Pro RUBIN (Loop Connection)  
Dementia-sensitive colour concept

## Accessories / Options

### NORA Sit-to-Stand Sling with Clips

| Canadian Orders Only |           | US Orders Only |             |
|----------------------|-----------|----------------|-------------|
| Description          | Part No.  | Description    | Part No.    |
| Small                | 921070053 | Small          | 921070053US |
| Medium               | 921070103 | Medium         | 921070103US |
| Large                | 921070203 | Large          | 921070203US |
| Extra Large          | 921070303 | Extra Large    | 921070303US |

### NORA Sit-to-Stand Sling with Loops

| Description                    | Part No.  |
|--------------------------------|-----------|
| Sling with Loops – Small       | 921071050 |
| Sling with Loops – Medium      | 921071100 |
| Sling with Loops – Large       | 921071200 |
| Sling with Loops – Extra Large | 921071300 |

### Spare Battery

| Description         | Part No.  |
|---------------------|-----------|
| Replacement battery | 921077400 |

### Wall Charger

| Description      | Part No.  |
|------------------|-----------|
| Charging station | 921077300 |

### Leg Separator

| Description   | Part No.  |
|---------------|-----------|
| Leg separator | 921077100 |

### HANDICARE

81 Romina Drive, Concord, ON, L4K 4Z9  
Toll Free: 877.304.5438 | Fax: 416.260.0256  
customerservice.canada@handicare.com  
www.handicare.ca

10888 Metro Court, St. Louis, MO 63043  
Toll Free: 866.891.6502 | Fax: 314.692.7858  
customerservice.northamerica@handicare.com  
www.handicareusa.com

Are you ready to make safe and effortless client transfers part of your everyday?

To learn more, call 1-888-637-8155 or visit [www.handicare.ca](http://www.handicare.ca) / [handicareusa.com](http://handicareusa.com)