



# PREMIUM PLUS Bathtub

The Perfect Bathing Solution

# Safety and Comfort All While Bathing

The PREMIUM PLUS Bathtub is designed with caregiver and resident safety in mind. Compatible with all bath lifts on the market with electronic height adjustment providing a low entry for independent transfers from a wheelchair if needed. Its ergonomic design provides the resident with a safe enjoyable bathing experience on its own or with the addition of optional features like the Air-Jet Whirlpool and Colour Light Therapy Systems. Made using Glass Reinforced Plastic (GRP) the PREMIUM PLUS is the perfect bathing solution with a modern forward thinking design.



## PREMIUM PLUS Bathtub Options




### **Air-Jet Whirlpool System**

Creates a pleasing effect on the skin, promoting a more relaxed tranquil bathing experience.




### **Colour Light Therapy System**

Transforms the bathing water into a vibrant colour spectrum creating a stimulating atmosphere.

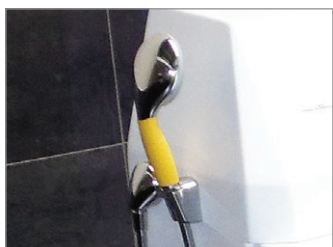

### **Disinfection System**

Offers increased protection against cross infection with approved in-house cleaning fluids.

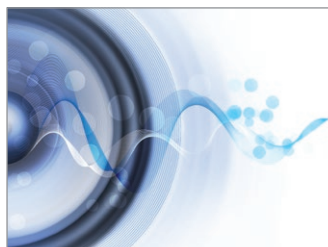

### **Music Sound System**

Underwater speakers generate gentle vibrations of music into the water aiding in relaxation.



### **Automatic Tub Filling**

Saves time for the caregiver by auto-filling the tub with warm water for immediate immersion.



### **Removable Tub Shortener**


For smaller residents providing more comfort and support.



# Caregiver and Resident Benefits



Digital operational panel includes three temperature displays: shower water, spout water and tub water. (optional)

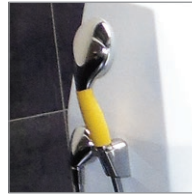

The disinfection system offers maximum protection against cross contamination. The disinfection hand shower and disinfectant container area are secured behind the control area away from the clients reach. (optional)



Made using Glass Reinforced Plastic (GRP) for increased durability.



Accessible from three sides as well as underneath. No need to be fixed to the floor.



Automatic water filling stop enables you to select and set water levels. (optional)



The lowered leg section, back profile and shaped armrests are ergonomically designed for resident comfort. The seat area also features a large hygiene opening for easy cleaning of the intimate area.



Chrome brass grab handles give the resident additional leverage and support.



Digital operational panel with visual symbols for easy identification and precision.

# PREMIUM PLUS Bathtub

## Specifications

### Safe Working Load (SWL)

Max patient weight: 220 kg (485 lb)

Patient and water: 300 kg (661 lb)

### Dimensions (L x W)

#### PREMIUM PLUS 175

Length: 2080 mm (81.9 in)

Width: 820 mm (32.3 in)

Height: 960 – 1360 mm (37.8 – 53.5 in)

#### PREMIUM PLUS 190

Length: 2230 mm (87.8 in)

Width: 820 mm (32.3 in)

Height: 960 – 1360 mm (37.8 – 53.5 in)

### Capacity

#### PREMIUM PLUS 175

Without patient: 190 kg (50.1 gal)

#### PREMIUM PLUS 190

Without patient: 220 kg (58.1 gal)

### Lifting Range

662 – 1062 mm (26.1 – 41.8 in)

### Water Drainage

45 L / minute (11.88 gal / minute)

### Filling

Volume, without client: appx. 220 L (58.1 gal)

Filling time: 42 L / minute (11.1 gal / minute)

### Power Supply

120V, 60 Hz

### Product Weight Without Packaging

175 kg (386 lb)

For additional technical data, please refer to the user manual.

## Part Numbers

**810142501:** PREMIUM PLUS 175

**810143501:** PREMIUM PLUS 190

## Accessories / Options

Description	PREMIUM PLUS 175 / 190
Tub shortener	810609200
Disinfection system	810328001
Disinfection blocking system	810328100
Temperature display for tub water	810146100
Automatic water stop	810131300
Double scalding protection function	810146100
Second hand shower	810328200
Music sound system	810131400
Air-jet whirlpool system	810131500
Colour light therapy system	810131600
Emergency function at power failure	810150857
Prestige safety grab bar	910303000



Are you ready to provide  
your residents with a unique  
relaxing bathing experience?

To learn more, call 1-888-637-8155 or  
visit [www.handicare.ca](http://www.handicare.ca) / [handicareusa.com](http://handicareusa.com)

### HANDICARE

485 Millway Avenue, Unit 2, Concord, ON, L4K 3V4

Toll Free: 877.304.5438 | Fax: 416.260.0256

[customerservice.canada@handicare.com](mailto:customerservice.canada@handicare.com)

[www.handicare.ca](http://www.handicare.ca)

10888 Metro Court, St. Louis, MO 63043

Toll Free: 866.891.6502 | Fax: 314.692.7858

[customerservice.northamerica@handicare.com](mailto:customerservice.northamerica@handicare.com)

[www.handicareusa.com](http://www.handicareusa.com)