

Span F350CT

Car Transfer Patient Lift



LIFTS



SURFACES



BEDS



SEATING



FALL PROTECTION



POSITIONERS



SKIN CARE



For users
Requiring
**MAXIMUM
ASSISTANCE**

- **The ideal solution** for simple transfers between car and chair, chair and bed, and from floor to chair.
- **Allows more dynamic positioning**, adjustment and control during transfer than with typical lift designs.
- **Portable and foldable** for easy transport and storage.
- **Safe, reliable, and easy operation:** 350 lb (160 kg) safe working load; simple charging, 3-year warranty.

SPAN
better science. better care.



F350CT Car Transfer Patient Lift

Created primarily for safe, controlled transfers into and out of a car—but with the flexibility for much more! Unique spreader bar design enables precise, controlled positioning with minimal motion.



Folds flat for easy transport and storage.




Foot pedal base width adjustment.



2-button pendant with LED indicator.



Underbed clearance allows use with most low beds.

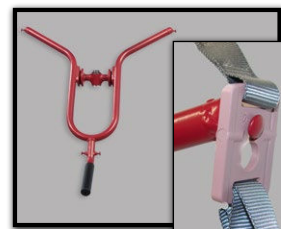

4-point spreader bar with "clip" sling attachment.

SPECIFICATIONS	
Weight Capacity	350 lb / 160 kg
External Base Width (closed)	27.6"
Internal Leg Width (opened)	32.0"
Mast Height	47.3"
Overall Base Length	41.0"
Minimum Lifting Height	22.2"
Maximum Lifting Height	64.7"
Net Weight	95 lb / 43 kg
Shipping Weight	108 lb / 49 kg

- **Locking rear casters**, foot-activated lock/unlock.
- **Floor extraction:** Can be used to lift patients from floor.
- **Quick-release battery** can be charged via charging cradle (optional), or via plug-in charger (included) for charging in place.
- **Overload protection** with visual/audible alert, halts lifting when capacity is exceeded.
- **Manual emergency lowering** for added safety.
- **Simple, fast assembly** using tools provided.
- **Certified to:** ISO 10535; IEC/EN/ANSI/AAMI/CAN/CSA-C22.2 60601-1, 60101-1-2, and 60601-1-11; SGS Mark for Safety.



- **LCD battery-level** indicator.
- **Audible low-battery** indicator.
- **Backlit** for easy viewing.
- **Redundant operation** via control unit; allows use if pendant lost/damaged.

Item #	Description
F350CT	Car Transfer Patient Lift
Warranty	3 years (frame and spreader bar); 2 years (electronics); 1 year (battery)

Part of Span's line of Safe Patient Handling and Mobility Solutions:



Full body and Sit-to-Stand floor lifts



Fixed and portable lithium-ion overhead lifts



Re-usable and disposable specialty slings