

# Span F400

Foldable / Portable Patient Lift



LIFTS



SURFACES



BEDS



SEATING



FALL PROTECTION



POSITIONERS



SKIN CARE



For users  
Requiring  
**MAXIMUM  
ASSISTANCE**




- **Folds flat simply and quickly** for transport and storage; fast, hassle-free assembly.
- **Safe and Reliable:** Redundant controls; overload protection; emergency lowering; 3-year frame warranty.
- **Easy operation:** Simple charging; audible/visual indicators; locking rear casters.
- **Powerful yet maneuverable:** 400 lb (182 kg) safe working load provides flexibility in a compact footprint.

**SPAN**  
better science. better care.



# F400 Foldable/Portable Patient Lift

*At home, on the road, or in-facility, the F400 is an outstanding choice for premium performance wherever space is at a premium.*

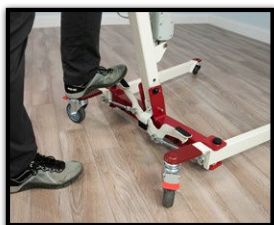


**Foot pedal base width adjustment.**



**2-button pendant with LED indicator.**



**Underbed clearance allows use with most low beds.**

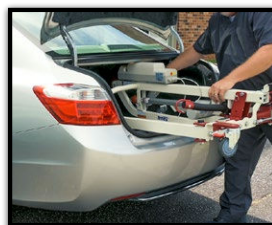

**Folds flat for easy transport and storage.**



**2-point spreader bar accommodates a variety of sling options.**

SPECIFICATIONS	
Weight Capacity	400 lb / 182 kg
Overall Base Length	41.0"
External Base Width (Closed)	27.6"
Internal Leg Width (Opened)	32.0"
Upright Mast Height	48.3"
Base Height	4.7"
Minimum Lifting Height with Spreader Bar	18.5"
Maximum Lifting Height with Spreader Bar	67.2"
Net Weight	90 lb / 41 kg
Shipping Weight	104 lb / 47 kg

- **Locking rear casters**, foot-activated lock/unlock.
- **Floor extraction:** Can be used to lift patients from floor.
- **Quick-release battery** can be charged via charging cradle (optional), or via plug-in charger (included) for charging in place.
- **Overload protection** with visual/audible alert, halts lifting when capacity is exceeded.
- **Manual emergency lowering** for added safety.
- **Simple, fast assembly** using tools provided.
- **Certified to:** ISO 10535; IEC/EN/ANSI/AAMI/CAN/CSA-C22.2 60601-1, 60101-1-2, and 60601-1-11; SGS Mark for Safety.

Item #	Description
F400	Foldable and Portable Patient Lift
Warranty	3 years (frame and spreader bar); 2 years (electronics); 1 year (battery)



- **LCD battery-level** indicator.
- **Audible low-battery** indicator.
- **Backlit** for easy viewing.
- **Redundant operation** via control unit; allows use if pendant lost/damaged.

## Part of Span's line of Safe Patient Handling and Mobility Solutions:



*Full body and Sit-to-Stand floor lifts*



*Fixed and portable lithium-ion overhead lifts*



*Re-usable and disposable specialty slings*

**SPAN**  
better science. better care.  
[www.spanamerica.com](http://www.spanamerica.com)