

Span F600B

Bariatric Full Body Patient Lift



*For users
Requiring
MAXIMUM
ASSISTANCE*

- **High Capacity:** 600 lb (273 kg) safe working load.
- **Care flexibility:** Size, maneuverability and range make it appropriate for both bariatric and non-bariatric users.
- **Safe, reliable and easy operation:** Simple charging, 3-year warranty.
- **Underbed clearance** allows use with most low beds.

SPAN
better science. better care.



LIFTS



SURFACES



BEDS



SEATING



FALL PROTECTION



POSITIONERS




SKIN CARE



F600B Bariatric Full Body Patient Lift

Heavy-duty performance in a manageable footprint. Delivers an unmatched package of high capacity combined with range, mobility, function, durability, ease of use and value not often found in bariatric models.

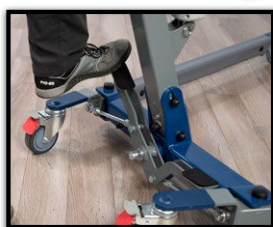


Foot pedal base width adjustment.




2-button pendant with LED indicator.



Underbed clearance allows use with most low beds.

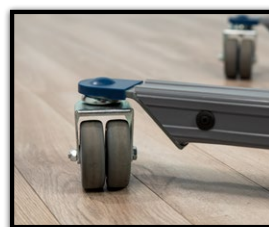

Dual front casters provide smooth support.

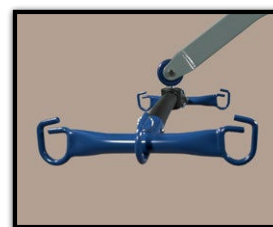

6-point spreader bar maximizes sling options.

SPECIFICATIONS	
Weight Capacity	600 lb / 273 kg
Overall Base Length	46.9"
External Base Width (Closed)	26.4"
Internal Leg Width (Opened)	39.3"
Upright Mast Height	51.3"
Base Height	4.7"
Minimum Lifting Height with Spreader Bar	21.5"
Maximum Lifting Height with Spreader Bar	68.6"
Net Weight	123 lb / 56 kg
Shipping Weight	139 lb / 63 kg

- **Locking rear casters**, foot-activated lock/unlock.
- **Floor extraction:** Can be used to lift patients from floor.
- **Quick-release battery** can be charged via charging cradle (optional), or via plug-in charger (included) for charging in place.
- **Overload protection** with visual/audible alert, halts lifting when capacity is exceeded.
- **Manual emergency lowering** for added safety.
- **Simple, fast assembly** using tools provided.
- **Certified to:** ISO 10535; IEC/EN/ANSI/AAMI/CAN/CSA-C22.2 60601-1, 60101-1-2, and 60601-1-11; SGS Mark for Safety.

Item #	Description
F600B	Bariatric Full Body Patient Lift
Warranty	3 years (frame and spreader bar); 2 years (electronics); 1 year (battery)



- **LCD battery-level** indicator.
- **Audible low-battery** indicator.
- **Backlit** for easy viewing.
- **Redundant operation** via control unit; allows use if pendant lost/damaged

Part of Span's line of Safe Patient Handling and Mobility Solutions:



Full body and Sit-to-Stand floor lifts



Fixed and portable lithium-ion overhead lifts



Re-usable and disposable specialty slings

SPAN
better science. better care.
www.spanamerica.com