


PressureGuard® Custom Care® Convertible LAL

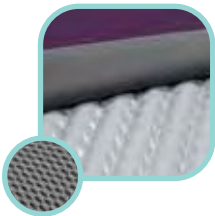

By its very design, the PressureGuard® Custom Care® Convertible LAL delivers improved patient outcomes for microclimate management and pressure redistribution with the added benefits associated with convertible surfaces.

The inner air delivery cover of the Custom Care Convertible LAL features a dedicated air supply and continuous airflow through its Air Diffusion Matrix™. This innovative design sweeps away moisture vapor before it can re-form as a liquid, helping to reduce maceration. Furthermore, with the powered control unit, the low air loss treatment helps to control excess moisture and removes excess heat at the microclimate level via moisture management.

Effective pressure management is possible through a number of patented design features, including our Geo-Matt® surface geometry and the gentle Heel Slope® of our dual heel protection.

¹ If the patient has a flap or graft that is less than 4 week post-op, total off-loading (proning) or air fluidized therapy is recommended. Avoid sliding patient (change resident to patient) on surface during repositioning, ADLs and transfers. Float or alternating pressure modes (not lateral rotation) are recommended when the patient is placed directly on the flap or graft. Source: *Wound Care Practice*. Sheffield, P. 1st edition, 2004, Chapter 17, P. 345. Other Span surfaces may be appropriate if used according to the guidelines.

² In situations where positioning on the pressure injury cannot be avoided (e.g. the individual has multiple injuries on multiple surfaces), limit the amount of time the individual is positioned on the injury. Source: *NPAUP/EPUAP Pressure Ulcer Prevention and Treatment, Clinical Practice Guideline*, p 66-67.

Therapy Provided

- ✓ Non-powered (Air) Reactive Pressure Redistribution
- With control unit added:*
- ✓ Microclimate Management (LAL)
- ✓ Alternating Pressure
- ✓ Lateral Rotation¹
- ✓ Powered Pressure Redistribution

QUICK FACTS

Indications:

- Excessive skin moisture due to sweating, decreased mobility, inability or refusal to be repositioned
- Macerated skin due to any of the above
- Increased skin or body temperature due to infection, sepsis, or other conditions
- A combined score ≤ 5 on Braden Moisture, Activity and Mobility scales
- Treatment of Stage 1-4 Pressure Injuries
- Prevention of Pressure Injuries for high-risk patients
- Some patients with stabilized flaps/grafts (see Algorithm, page 57)¹

Key Advantages:

- Effective pressure transfer through Geo-Matt® surface geometry and dual heel protection
- Reduced maceration through continuous airflow via Air Diffusion Matrix™ layer
- More effective powered therapies through the Star Chamber™ design (see pages 66-67)
- Easy-to-use control unit with quick-connect air lines and push-button therapy selection

Item #	Description
CL753629	75"L x 36"W x 7"H
CL803629	80"L x 36"W x 7"H
CL843629	84"L x 36"W x 7"H
74750	80"L x 42"W x 7"H
49907	84"L x 42"W x 7"H
* CL85V29	Custom Care Convertible LAL for recessed deck frames (see page 17)

* Allow up to 10 working days lead time for this version.

Control Unit:	8410	Add-on control unit for Custom Care Convertible LAL
Medicare Approved:	E0373 (Group II) (without control unit) E0277 (Group II) (with control unit)	
Weight limit:	Supports up to 500 lbs. (non-powered & alternating modes) Supports up to 350 lbs. (lateral rotation mode)	
Weight:	Mattress: approximately 34 lbs. Control unit: 6 lbs.	
Warranty:	Mattress: 5 years, non-prorated Control unit: 2 years, non-prorated	
UL listed:	Yes	
Flammability:	All models are certified for conformance to NFPA 101 (Life Safety Code) – ASTM E1590; Cal Tech Bulletin #117 and 16 CFR Parts 1632 and 1633 by an independent testing organization.	

