


Geo-Mattress® UltraMax™ for Stretchers

Category:
Premium Pressure Redistribution
Support Surface





3-Dimensional Zoning™ design provides anatomically-calibrated levels of support and immersion throughout.

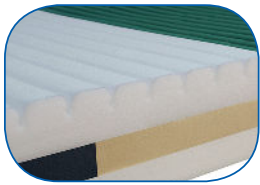

Extra-soft heel section facilitates immersion to provide added heel protection.

Design Advantages Inside...

Engineered, multi-layer support core

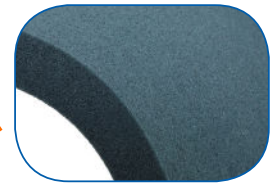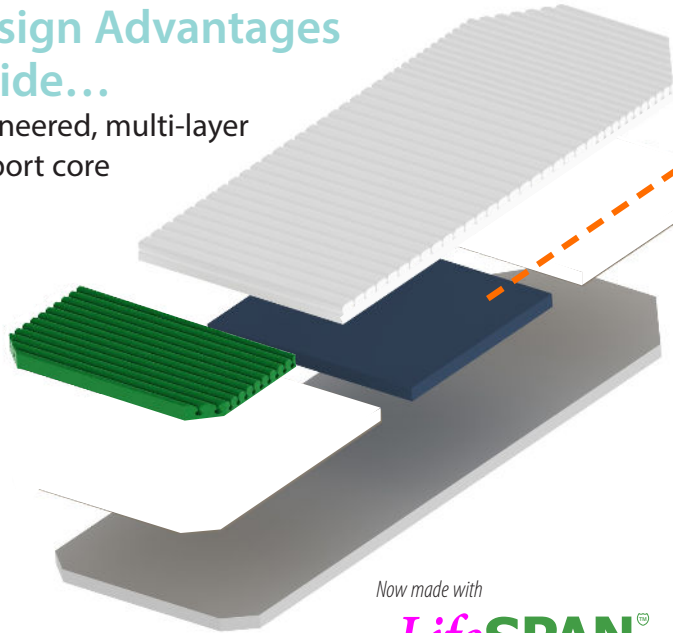

UHP 500™ Ultra High Performance foam in high-fatigue seat section; 55% higher support factor yields superior resiliency and recovery.



Geo-Matt® non-shearing surface creates individually articulating cells; Ring-of-Air® channels dissipate heat and moisture.

Now made with



bleach-resistant fabric cover
(see page 15 for details)

Specifications

[All models listed are 5" high to facilitate maximum immersion]

Item #	Dimensions, overall (L)ength and (W)idth	(H)ead-end taper length	(F)oot-end taper length
47318	75"L x 26"W	7"	0"
49297	75"L x 26"W	13"	4"
97107	76"L x 26"W	10"	4"
10635	76"L x 26"W	12"	5"
2879529	76"L x 27"W	11"	5"
46234	76"L x 29"W	10"	8"
79012	76"L x 30"W	13"	4"
98969	78"L x 27"W	12"	12"
81921	2-pc "Gynie"	n/a	n/a

Other dimensions also available upon request.


Supports: up to 700 lbs.

Warranty: 3-year, non-prorated

Mattress Weight: approx. 20 lbs.

All components: 100% latex-free and 100% PBDE-free

Flammability: All models are certified for conformance to NFPA 101 (Life Safety Code) – ASTM E1590; Cal Tech Bulletin #117 and 16 CFR Parts 1632 and 1633 by an independent testing organization.

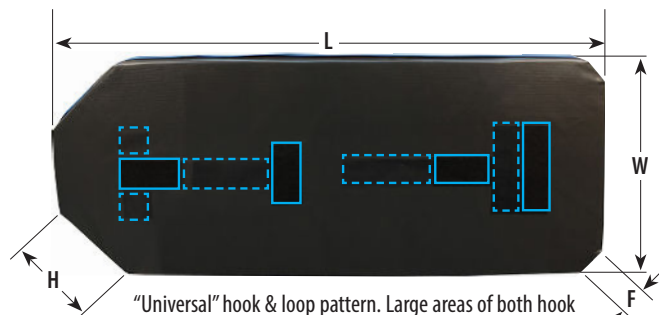

"Universal" hook & loop pattern. Large areas of both hook and loop on bottom of mattress help in anchoring it to a wide variety of stretcher models.

WARNING: Cancer and Reproductive Harm – www.P65Warnings.ca.gov



SURFACES