Geo-Mattress with Wings®





Item #	Description
W7535-29	75"L x 35"W x 6"H center, 8"H side
W7835-29	78"L x 35"W x 6"H center, 8"H side
W8435-29	84"L x 35"W x 6"H center, 8"H side









angled bolster, topped with proprietary Geo-Matt® segmented surface.

Stable, high-density

Geo-Mattress with Wings®

Safety — without compromise. Helps protect skin integrity while providing gentle prompting away from the edge of the bed.

- 5-year, non-prorated warranty
- · Weight Limit: 400 lbs.

Flammability Protection: Certified conformance to NFPA 101 (Life Safety Code); Cal. Tech Bulletin #117; 16 CFR 1632 and 1633

Primary Use:

• Prevention of Skin Breakdown


Patient Indications:

- · Patients at fall risk due to independent exits from bed
- · General patient population
- Moderate to high risk for skin breakdown
- · Decreased mobility
- Repositions frequently by self or caregiver


Key Advantages:

- Edge-to-edge antishearing surface
- Extra-firm edge bolster
- · Shearing protection at transfer opening
- Heel Slope[®]

Packed with value-added features:

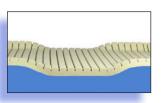

Lowered, padded openings make transfers easier and safer. Proprietary edge-to-edge Geo-Matt segmented surface provides shearing protection.

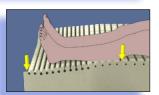


Exclusive Heel Slope™.

Subtle taper redistributes load to pressure-tolerant lower legs. Delivers heel pressures 27% lower than flat designs.



Fluid-resistant, wipe-clean cover with ultra-low moisture vapor transmission rate (MVTR). Bacteriostatically-treated top, durable tri-laminate bottom. Hidden zipper design for easier, more thorough cleaning and superior infection control.

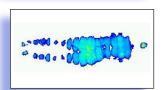

Effective pressure redistribution across the surface ensures low average of peak pressure at vulnerable sacrum and heels.

